

SPX-29

SCALERS, PERIODONTOLOGY

Efficiency. Precision
Care — gentle removal,
superior control.

Surgiplex. The Difference You Hold



www.surgiplexsrl.com

Scalers, supragingival

Scaler, supragingivaler Einsatz

Limpiadores de sarro, supragingival

Instruments à détartrer, supragingival

Ablatori di tartaro, supragingival


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología


Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

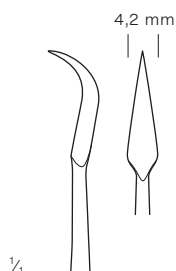


SP29-116-01

Fig. 1

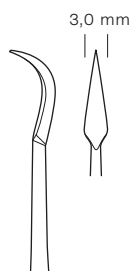


SP29-116-02

Fig. 2

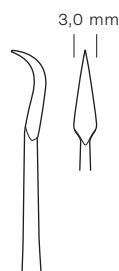


SP29-116-03

Fig. 3

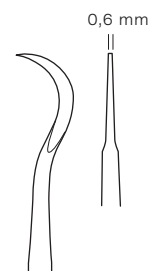


SP29-116-32

Fig. 32

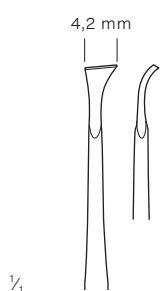

SP29-116-33

Fig. 33

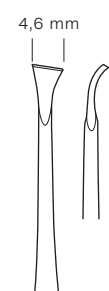

SP29-116-34

Fig. 34


**SP29-116-01 -
SP29-116-34**

17,5 cm / 6 7/8"


solid handle, 6 mm




**SP29-117-33 -
 SP29-121-07**
 17,5 cm / 6 7/8"
 solid handle, 6 mm

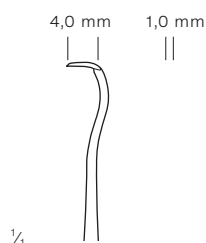


**Jaquette
 SP29-117-33**
 Fig. 33

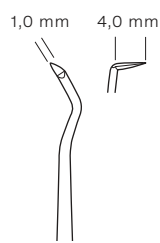


**Jaquette
 SP29-117-34**
 Fig. 34

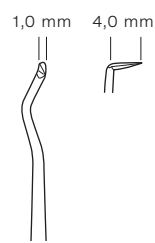


**Jaquette
 SP29-117-35**
 Fig. 35

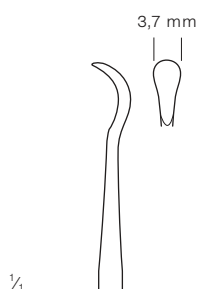


**Cumine
 SP29-118-04**
 Fig. 33

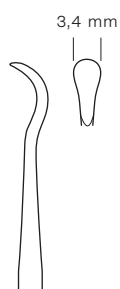




**Mitchell
 SP29-119-04**
 Fig. 34

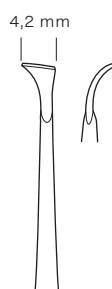

SP29-120-33
 Fig. 33/34

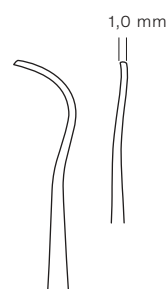

**Younger-Good
 SP29-121-07**
 Fig. 33/34

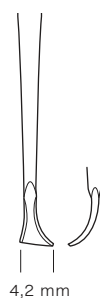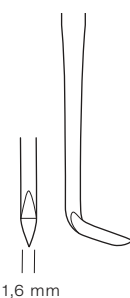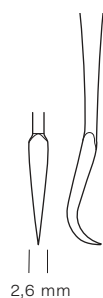

Scalers, supragingival

Scaler, supragingivaler Einsatz

Limpiadores de sarro, supragingival

Instruments à détartrer, supragingival

Ablatori di tartaro, supragingival

Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia




Towner
SP29-122-05

Fig. 5/33



Incisors, premolars,
universal, also approximal

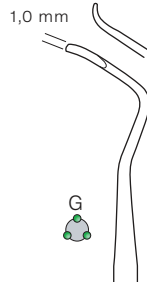

McCall
SP29-122-13

Fig. 13S/14S



Lateral teeth, mesial/distal
Particularly useful for working in
interdental spaces in molar region




McCall
SP29-122-14

Fig. 13SM/14



Lateral teeth; mesial/distal,
like 13S/14S, with rounded tip

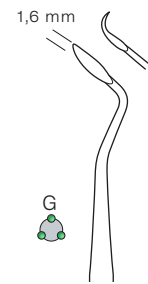

Taylor
SP29-122-15

Fig. T 2/3



Incisors, premolars, molars,
universal
Particularly useful for working
in interdental spaces in molar
region

SP29-122-05 -
SP29-122-15

17,5 cm / 6 7/8"

hollow handle, 8 mm

GR = grey

Y = yellow

R = red


G = green

V = violet

W = white

B = blue


N = black

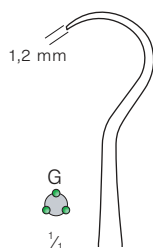

Towner
SP29-122-16
Fig. 15/30




For removing stubborn
tartar/scale

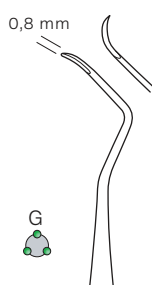

SP29-122-17
Fig. U 135



for hard-to-access
approximal surfaces

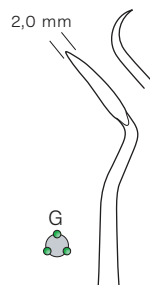


SP29-122-18
Fig. CI 2/3




Incisors, premolars, molars,
universal



SP29-122-19
Fig. M 23

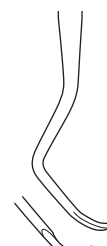

Incisors, premolars,
mesial/distal

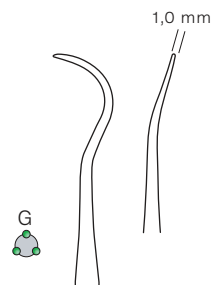

Mc Call
SP29-123-18
Fig. 17/18



universal

SP29-122-16 -
SP29-123-18

17,5 cm / 6 7/8"

hollow handle, 8 mm

Scalers, supragingival

Scaler, supragingivaler Einsatz

Limpiadores de sarro, supragingival

Instruments à détartrer, supragingival

Ablatori di tartaro, supragingival


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología


Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

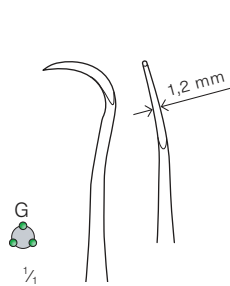


Hygienist
SP29-124-06

Fig. H6/H7

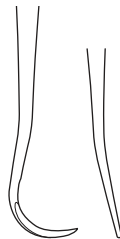

Incisivos, premolares,
mesial/distal
Esp. indicado para los
espacios interdentes

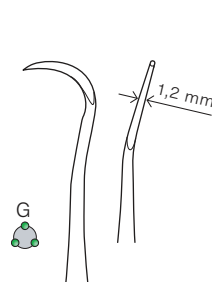

Hygienist
SP29-124-10

Fig. H6/H7



Incisivos, premolares,
mesial/distal
Esp. indicado para los
espacios interdentes



SP29-125-00

Fig. 204



Molars, universal



SP29-125-01

Fig. 204 S



Molars, universal



SP29-125-02

Fig. 204 SD



Incisors, premolars,
mesial/distal

SP29-124-06 -

SP29-125-02

17,5 cm / 6 3/4"

hollow handle, 8 mm

GR = grey

Y = yellow

R = red


G = green

V = violet

W = white

B = blue

N = black

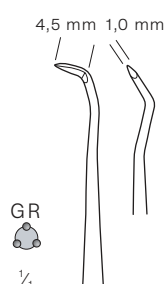

Zürich
SP29-126-17



For cleaning cavities
following apicœctomy



SP29-127-17
Fig. CC 6



SP29-128-17
Fig. T 2/3



Incisors, premolars,
molars, universal

**SP29-126-17-
SP29-128-17**
17,5 cm / 6 3/4"
hollow handle, 8 mm

Scalers, supragingival

Scaler, supragingivaler Einsatz

Limpiadores de sarro, supragingival

Instruments à détartrer, supragingival

Ablatori di tartaro, supragingival


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología


Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

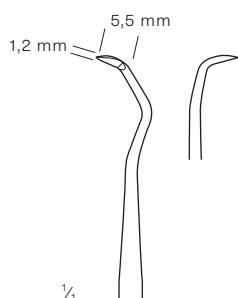


**Jaquette
SP29-129-32**

Fig. 31/32

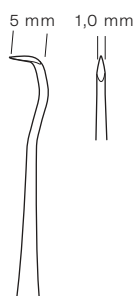

**Jaquette
SP29-129-33**

Fig. 30/33

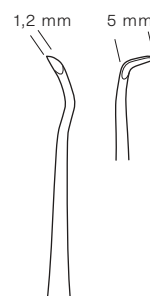

**Jaquette
SP29-129-35**

Fig. 34/35



SP29-129-32 -


SP29-129-35

17,5 cm / 6 7/8"


hollow handle, 8 mm




**SP29-130-31 -
SP29-130-36**
17,5 cm / 6 7/8"
hollow handle, 8 mm

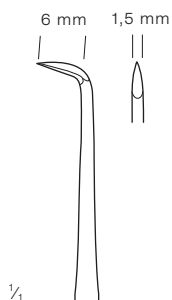


**Goldman-Fox
SP29-130-31**
Fig. 1

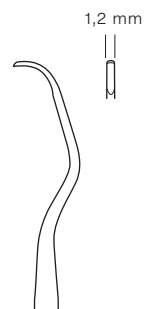


**Goldman-Fox
SP29-130-32**
Fig. 2

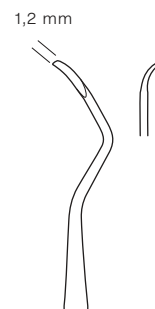


**Goldman-Fox
SP29-130-33**
Fig. 3

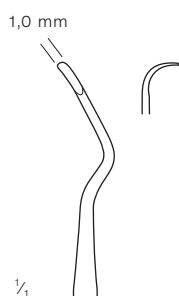

**Goldman-Fox
SP29-130-34**
Fig. 4

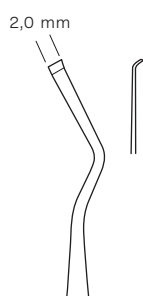

**Goldman-Fox
SP29-130-35**
Fig. 5

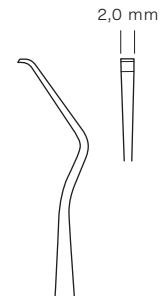

**Goldman-Fox
SP29-130-36**
Fig. 6



Scalers, supragingival, bone files

Scaler, supragingivaler Einsatz, Knochenfeilen
 Limpiadores de sarro, supragingival, limas de hueso
 Instruments à détartrer, supragingival, râpes à os
 Ablatori di tartaro, supragingival, lime e raspe per osso

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
 Limpiadores de sarro, parodontología
 Instruments à détartrer, parodontologie
 Ablatori di tartaro, parodontologia




**Hirschfeld
SP29-131-03**
Fig. 3/7




**Hirschfeld
SP29-131-05**
Fig. 5/11



File
 For removal of granulosomatous tissue from
 supra- and infrabone pockets

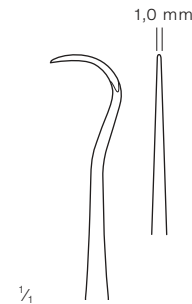

**Mercitan
SP29-132-05**
Fig. SH 5

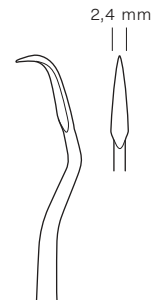

**Ward-Kirkland
SP29-132-10**
Fig. 1



**Crane-Kaplan
SP29-132-20**
Fig. 2



**Crane-Kaplan
SP29-132-30**
Fig. 3




**White
SP29-133-12**
Fig. 12B

**SP29-131-03 -
SP29-131-05**
17,5 cm / 6 7/8"
hollow handle, 8 mm


**SP29-132-05 -
SP29-132-30**
17,5 cm / 6 7/8"
solid handle, 6 mm




SP29-134-01 -
SP29-134-08
17,5 cm / 6 7/8"
hollow handle, 8 mm

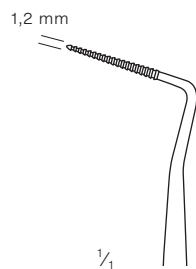



SP29-134-01

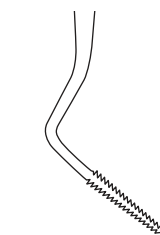


modified

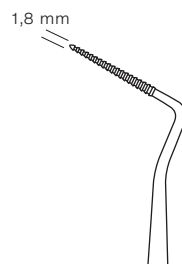

SP29-134-02

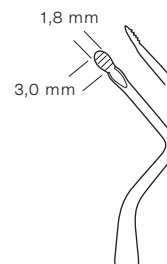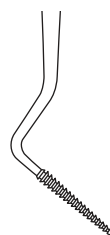


Orban
SP29-134-06
Fig. 12/13

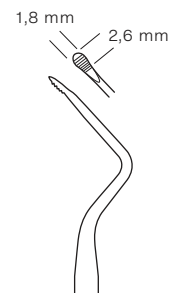

SP29-134-08
Fig. 12/12




For removal of granulosomatous tissue from supra- and infrabone pockets

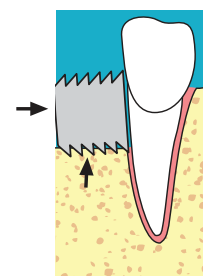

Periodontal file
 Parodontalfeile
 Lima parodontal
 Lime parodontale
 Lima periodontale


Scalers, periodontology
 Reinigungsinstrumente, Parodontologie
 Limpiadores de sarro, parodontología
 Instruments à détartrer, parodontologie
 Ablatori di tartaro, parodontologia

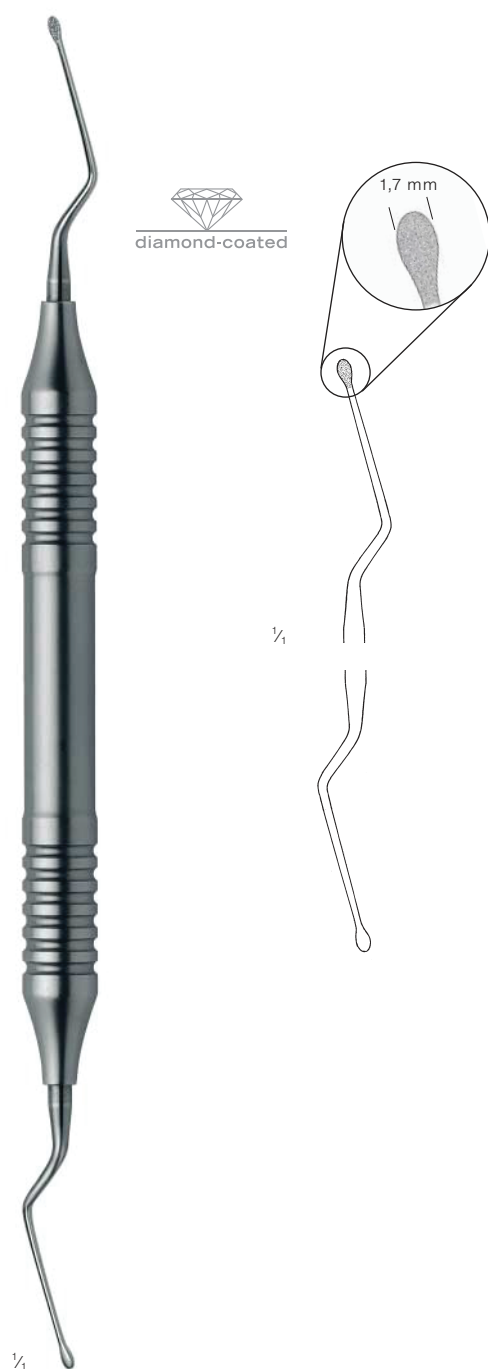


Roncati
SP29-135-00
 17,5 cm / 6 7/8"
 hollow handle, 10 mm



Roncati
SP29-135-01
 17,5 cm / 6 7/8"
 hollow handle, 10 mm

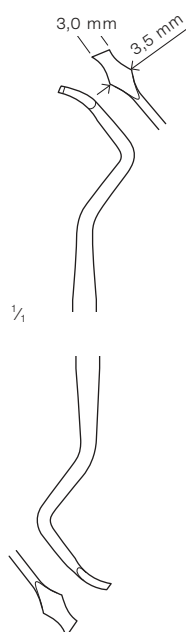

Mini chisel

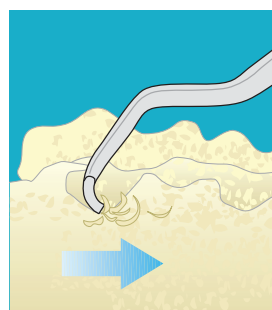

Rhodes
SP29-136-01


Fig. 36/37

17,5 cm / 6 7/8"

hollow handle, 8 mm


Periodontal chisels
 Parodontalmeißel
 Cinceles de periodoncia
 Burins de parodontie
 Sgorbie per parodontologia

Scalers, periodontology
 Reinigungsinstrumente, Parodontologie
 Limpiadores de sarro, parodontología
 Instruments à détartrer, parodontologie
 Ablatori di tartaro, parodontologia

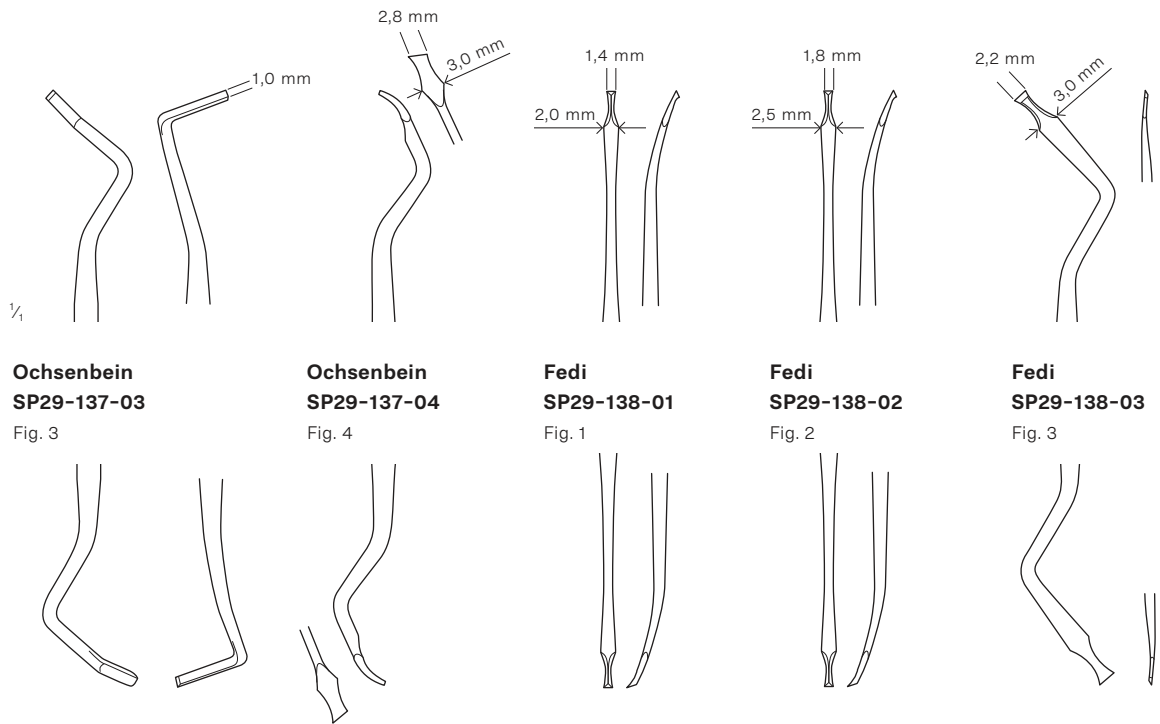

For flap operation. Also for minor bone corrections


**SP29-137-03 -
 SP29-138-03**
 17,5 cm / 6 7/8"
 hollow handle, 8 mm

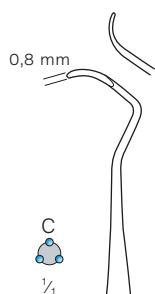

Columbia
SP29-139-13
 Fig. 13/14




Incisors and lateral teeth,
 universal

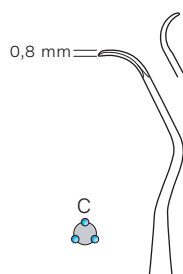

Columbia
SP29-139-20
 Fig. 2R/2L




Incisors, universal

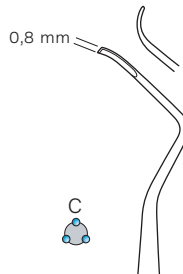


Columbia
SP29-139-40
 Fig. 4R/4L




Lateral teeth, universal

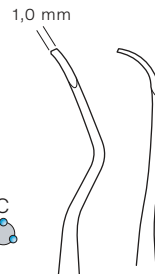

SP29-140-02
 Fig. GX4

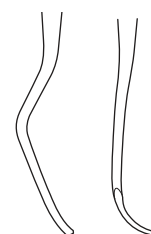

Incisors and lateral teeth,
 universal

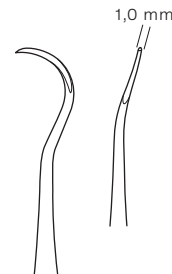

Mc Call
SP29-141-02
 Fig. 17S/18S



universal
 bifurcate root surfaces

SP29-139-13 -
SP29-141-02
 17,5 cm / 6 7/8"
 Universal curette
 hollow handle, 8 mm

GR grey Y yellow R red G green V violet W white B blue C light-blue

Curettes, subgingival

Küretten, subgingivaler Einsatz

Curetas, subgingival

Curettes, subgingival

Curette, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

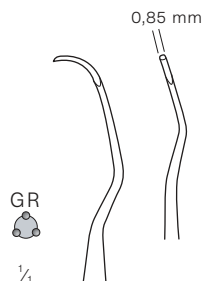


Langer	Langer
SP29-142-02	SP29-142-04
SP29-143-02	SP29-143-04

Fig. 1/2



Mandible, lateral teeth,
mesial / distal

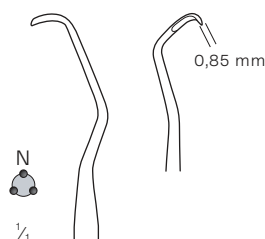

Langer	Langer
SP29-142-04	SP29-142-06
SP29-143-04	SP29-143-06

Fig. 3/4



Maxilla, lateral teeth,
mesial / distal



Langer	Langer
SP29-142-06	SP29-142-08
SP29-143-06	SP29-143-08

Fig. 5/6



Maxilla and mandible,
incisors, universal

Langer
SP29-142-02 -
SP29-143-06

17,5 cm / 6 7/8"

Universal curette

hollow handle, 8 mm

GR = grey

Y = yellow

R = red

G = green

V = violet

W = white

B = blue

N = black




Roncati Scaler
SP29-144-00
 17,5 cm / 6 7/8"
 Scaler
 hollow handle, 10 mm



Roncati
SP29-144-01
 17,5 cm / 6 7/8"
 Curette, micro, straight
 hollow handle, 10 mm



Roncati
SP29-144-02
 17,5 cm / 6 7/8"
 Curette, micro, curved
 hollow handle, 10 mm

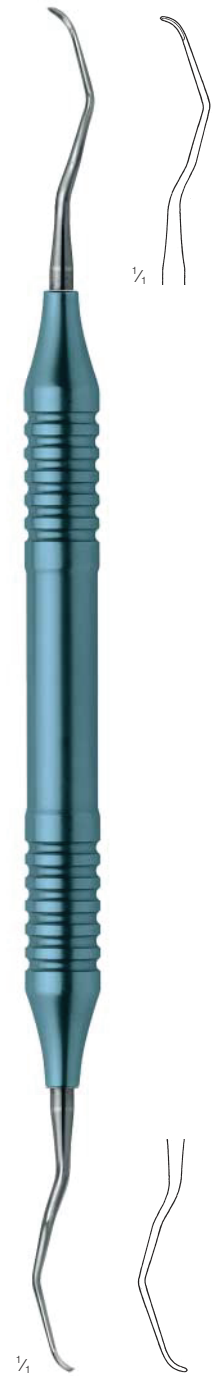

Roncati
SP29-145-00
 17,5 cm / 6 7/8"
 Curette, micro, curved
 hollow handle, 10 mm

Ti

Scaling- Instruments

Scaling Instruments
Instruments de Détartrage
Scaling-Instrumente

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia

Why scaling instruments should always be sharp:


1. Sharp instruments allow precision work with little application of force.
2. Blunt instruments polish the dental stone, bruise the soft tissue and require more application of force.

What you need for instrument sharpening:

1. Sharpening stones (rough, fine, and round for burr removal)
2. Sharpening oil (prevents excessive heat build-up; preserves the sharpening capacity of the sharpening stone)
3. Test stick (made of acrylic, for sharpness test)
4. Good light (sharp edges do not reflect light)
5. Clamping device (e.g. vise)

What to check after sharpening:

1. Roundness of back surfaces
2. Continuous sharpness of cutting edges
3. Perfect shape of working end (e.g. check against catalog illustration)
4. After several sharpening cycles, verify that the instrument is still strong enough for application. Replace instrument if necessary.

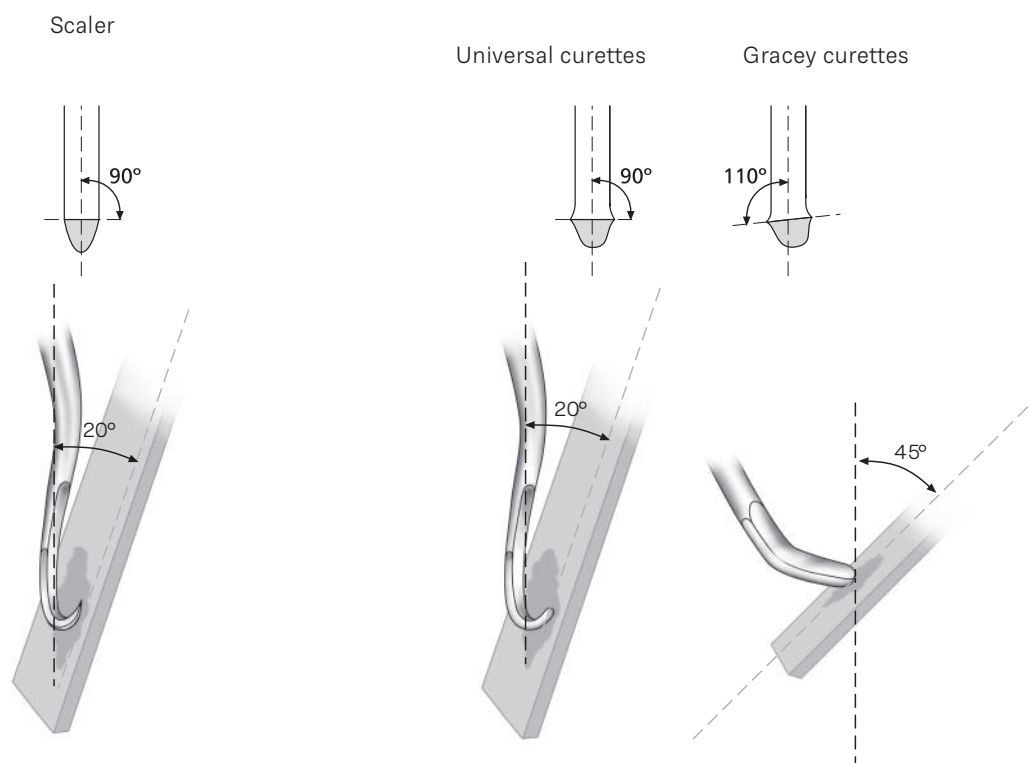

Color coding of scalers and Gracey curettes

1. All scalers are color-coded green – 3 points distributed across the diameter. At the same time, the color code marks the finer working end.
2. The Gracey Standard, Gracey Rigid and Gracey Mini curettes also have the color code on the side of the finer working end.
Example: Gracey 1/2 = color code applied on side of working end no. 1.
3. For even better differentiation between Gracey Standard and Gracey Mini on the one hand and Gracey Rigid on the other hand, Gracey Rigid curettes are coded with six color points.
4. Color-coded Gracey Standard curettes are also available with a 10 mm handle.


left: color
(lower figure,
e.g. insert 3)



right: no color
(higher figure,
e.g. insert 4)



	Ø		
Gracey Mini	8 mm	x	
Gracey Standard	8 mm	x	
Gracey Big	10 mm	x	
Gracey Rigid	8 mm		x

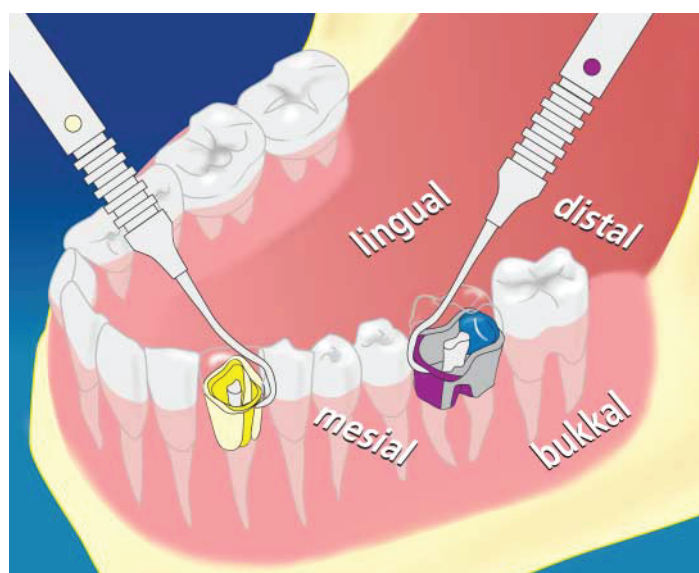

Fields of application for Gracey curettes

-  1-2,3-4, 5-6 Incisors
-  7-8, 9-10 Premolars and molars (buccal/lingual)
-  11-12, 15-16 Premolars and molars, mesial Furcation: mesial portion of the distal root
-  13-14, 17-18 Premolars and molars, distal Furcation: distal portion of the mesial root

Gracey mini curettes, subgingival

Gracey Mini-Küretten, subgingivaler Einsatz

Gracey mini curetas, subgingival

Gracey mini curettes, subgingival

Gracey mini curette, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

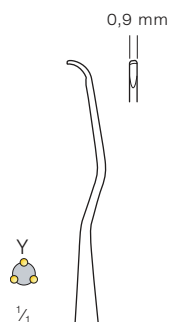


SP29-146-02

Fig. 1/2



Incisors, universal,
first shaft longer, working
end shorter

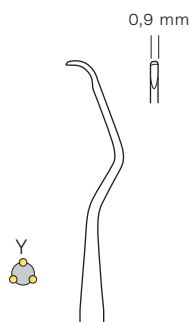


SP29-146-04

Fig. 3/4



Incisors, premolars,
universal,
first shaft longer,
working end shorter

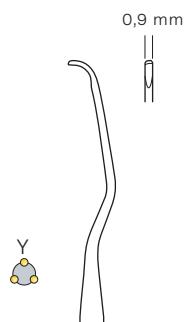


SP29-146-06

Fig. 5/6




Incisors, premolars,
universal,
first shaft longer,
working end shorter

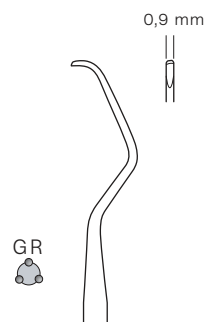

SP29-146-08

Fig. 7/8

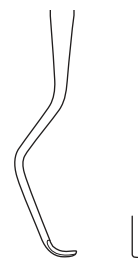

Premolars, molars,
lingual/buccal,
first shaft longer, working
end shorter

1/1

Gracey Mini

SP29-146-02 -

SP29-146-18

17,5 cm / 6 7/8"

hollow handle, 8 mm

GR = grey

Y = yellow

R = red


G = green

V = violet

W = white


B = blue

N = black

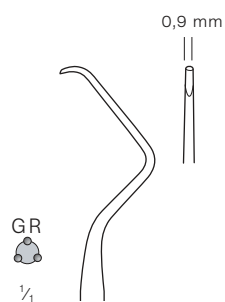


SP29-146-10

Fig. 9/10

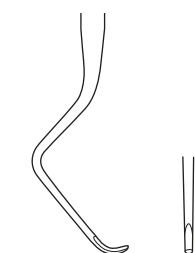


Premolars, molars,
lingual/buccal,
first shaft longer,
working end shorter

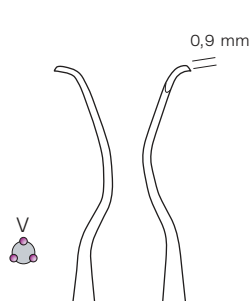


SP29-146-12

Fig. 11/12

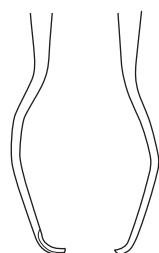


Premolars, molars,
mesial, mesiolingual,
mesiobuccal,
first shaft longer,
working end shorter

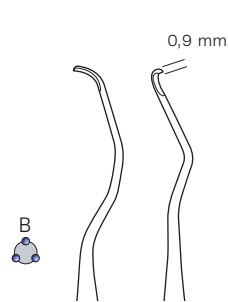


SP29-146-14

Fig. 13/14

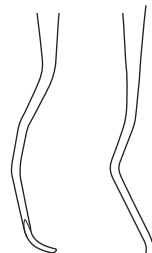


Premolars, molars, distal,
first shaft longer,
working end shorter

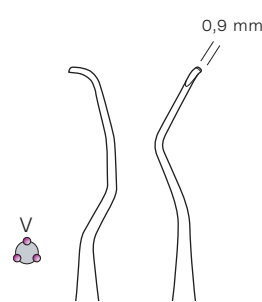


SP29-146-16

Fig. 15/16

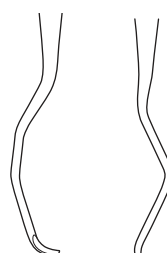


Premolars, molars, mesial,
first shaft longer,
working end shorter

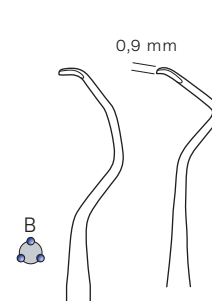

SP29-146-18

Fig. 17/18

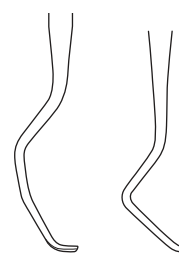

Premolars, molars, distal,
first shaft longer,
working end shorter

Gracey standard curettes, subgingival

Gracey Standard-Küretten, subgingivaler Einsatz

Gracey curetas estándar, subgingival

Gracey curettes standard, subgingival

Gracey curette standard, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

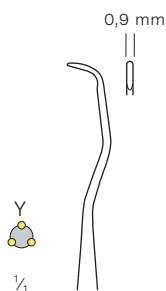


SP29-147-02

Fig. 1/2



Incisors, universal

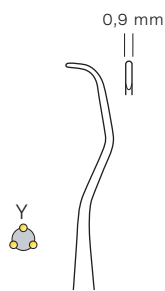


SP29-147-04

Fig. 3/4



Incisors, premolars,
universal

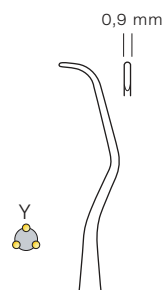


SP29-147-06

Fig. 5/6



Incisors, premolars,
universal

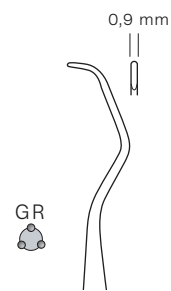

SP29-147-08

Fig. 7/8



Premolars, molars,
lingual / buccal

Gracey Standard


SP29-147-02 -

SP29-147-18


17,5 cm / 6 7/8"

hollow handle, 8 mm




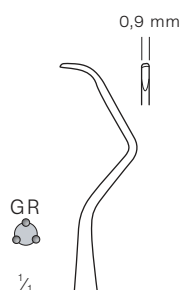

SP29-147-10
Fig. 9/10

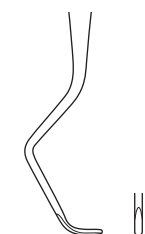


Premolars, molars,
lingual / buccal

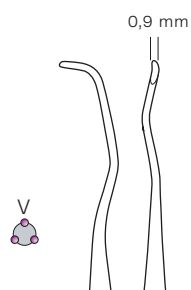

SP29-147-12
Fig. 11/12




Premolars, molars, mesial,
mesiolingual, mesiobuccal

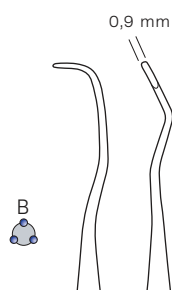

SP29-147-14
Fig. 13/14




Premolars, molars, distal

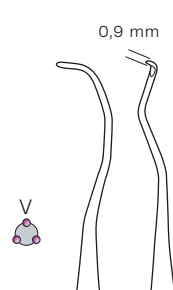


SP29-147-16
Fig. 15/16



Premolars, molars, mesial

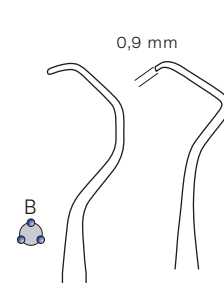

SP29-147-18
Fig. 17/18

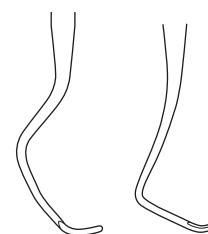

Premolars, molars, distal

Gracey Rigid curettes, subgingival

Gracey Rigid-Küretten, subgingivaler Einsatz

Gracey Rigid curetas, subgingival

Gracey Rigid curettes, subgingival

Gracey Rigid curette, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

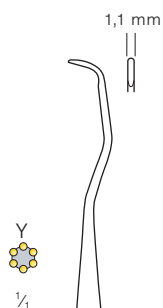


SP29-148-02

Fig. 1/2



Incisors, universal,
extra rigid, for stubborn
dental plaque

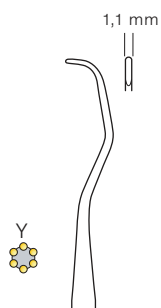


SP29-148-04

Fig. 3/4



Incisors, premolars,
universal,
extra rigid, for stubborn
dental plaque

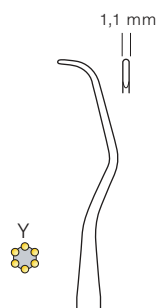


SP29-148-06

Fig. 5/6



Incisors, premolars,
universal,
extra rigid, for stubborn
dental plaque

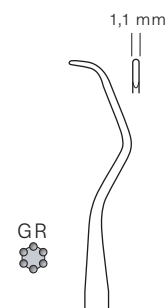

SP29-148-08

Fig. 7/8




Premolars, molars, lingual,
buccal,
extra rigid, for stubborn
dental plaque

Gracey Rigid SP29-148-02 - SP29-148-18

17,5 cm / 6 7/8"


hollow handle, 8 mm



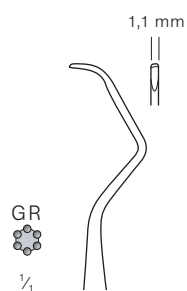


SP29-148-10

Fig. 9/10

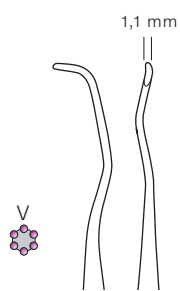


SP29-148-12

Fig. 11/12

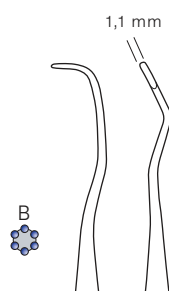


SP29-148-14

Fig. 13/14

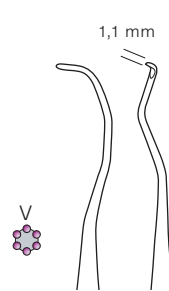


SP29-148-16

Fig. 15/16

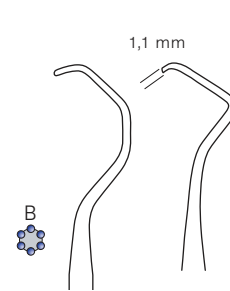

SP29-148-18

Fig. 17/18

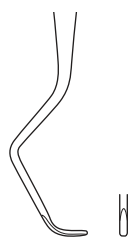

Premolars, molars, lingual,
 buccal,
 extra rigid, for stubborn
 dental plaque




Premolars, molars, mesial,
 mesiolingual, mesiobuccal,
 extra rigid, for stubborn
 dental plaque



Premolars, molars,
 distal
 extra rigid, for stubborn
 dental plaque



Premolars, molars, mesial,
 extra rigid, for stubborn
 dental plaque

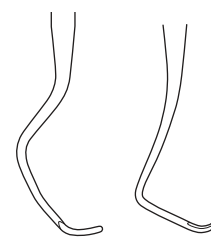


Premolars, molars,
 distal
 extra rigid, for stubborn
 dental plaque

Gracey curettes with 10 mm handle

Gracey Küretten mit 10 mm-Griff
 Gracey curetas para mangos 10 mm
 Gracey curettes avec manches à 10 mm
 Gracey curette con manici 10 mm

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
 Limpiadores de sarro, parodontología
 Instruments à détartrer, parodontologie
 Ablatori di tartaro, parodontologia

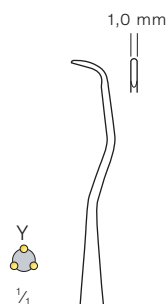


SP29-149-02

Fig. 1/2



Incisors, universal

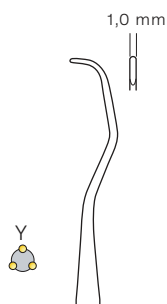


SP29-149-04

Fig. 3/4



Incisors, premolars, universal

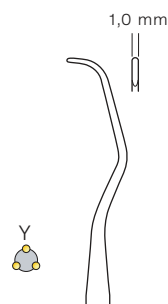


SP29-149-06

Fig. 5/6



Incisors, premolars, universal

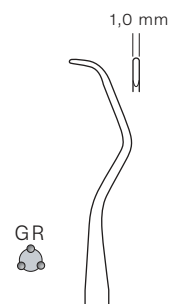

SP29-149-08

Fig. 7/8



Premolars, molars, lingual, buccal


Big Gracey

**SP29-149-02 -
 SP29-149-18**

17,5 cm / 6 7/8"


hollow handle, 10 mm



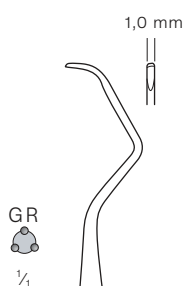


SP29-149-10

Fig. 9/10

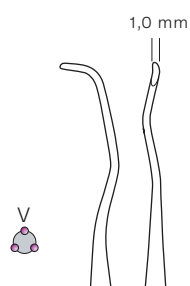


SP29-149-12

Fig. 11/12

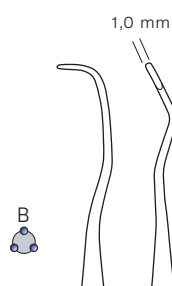


SP29-149-14

Fig. 13/14

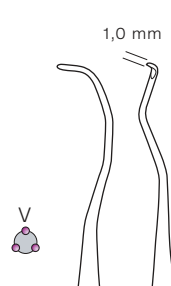


SP29-149-16

Fig. 15/16

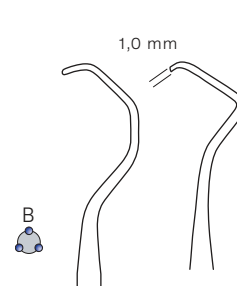

SP29-149-18

Fig. 17/18

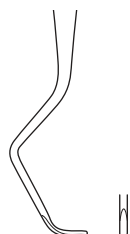

Premolars, molars, lingual,
buccal




Premolars, molars,
mesial, mesiolingual,
mesiobuccal,



Premolars, molars,
distal



Premolars, molars, mesial

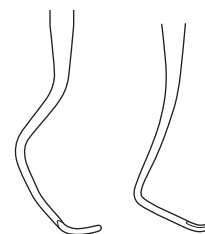

Premolars, molars,
distal

Curettes, subgingival

Küretten, subgingivaler Einsatz

Curetas, subgingival

Curettes, subgingival

Curette, subgingivale

Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia



SP29-150-01

Fig. 1



Furcation curette
0,9 mm
buccal / lingual



SP29-150-02

Fig. 2



Furcation curette
3 mm
buccal / lingual



SP29-150-03

Fig. 3



Furcation curette
0,9 mm
mesial / distal



SP29-150-04

Fig. 4




Furcation curette
1,3 mm
mesial / distal

**SP29-150-01 -
SP29-150-04**

17,5 cm / 6 7/8"

hollow handle, 8 mm

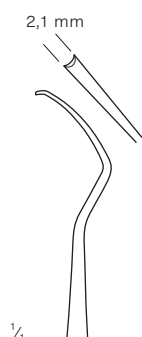


SP29-150-05

Fig. 5



Furcation curette
 0,9 mm
 labial / lingual

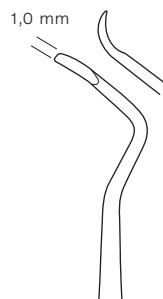

Prichard
SP29-151-01

Fig. 1/2



Bone curette,
 extra long, robust
 construction, double
 cutting

Gracey mini curettes, subgingival

Gracey Mini-Küretten, subgingivaler Einsatz

Gracey mini curetas, subgingival

Gracey mini curettes, subgingival

Gracey mini curette, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

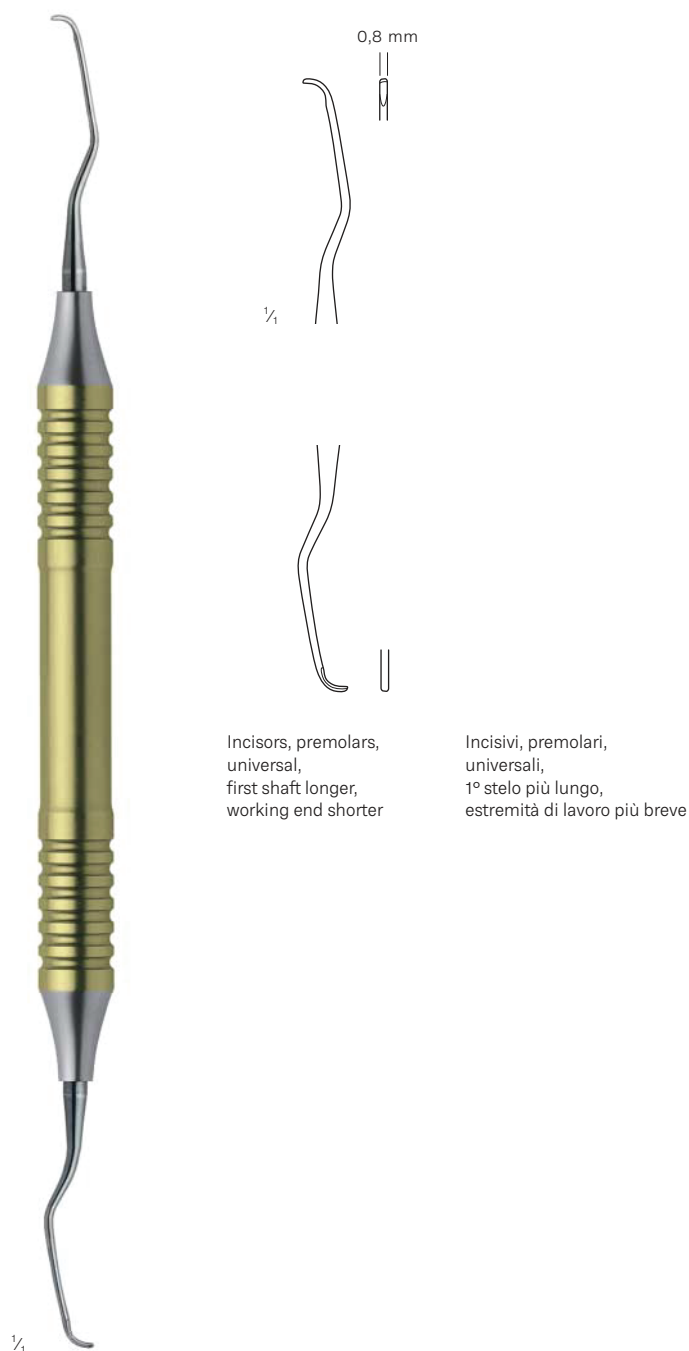


SP29-152-06

Ti

Fig. 5/6

17,5 cm / 6 7/8"

hollow handle, 10 mm

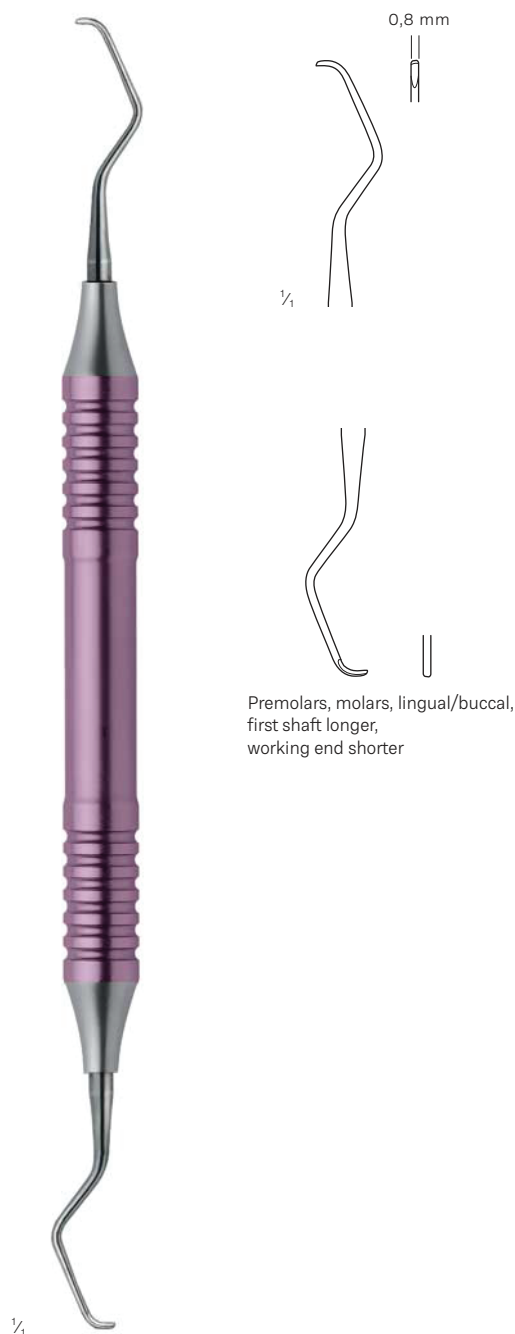



SP29-152-08

Ti

Fig. 7/8

17,5 cm / 6 7/8"

hollow handle, 10 mm

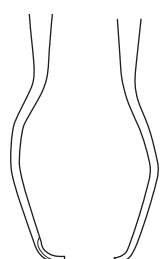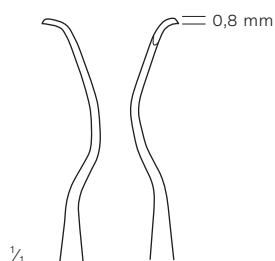

Premolars, molars,
mesial, mesiolingual,
mesiobuccal,
first shaft longer, working
end shorter

Premolari, molari,
mesiali, mesiolinguali,
mesiobuccali,
1° stelo più lungo, estremità
di lavoro più breve

1/1


SP29-152-12

Fig. 11/12

17,5 cm / 6 7/8"

hollow handle, 10 mm

Ti

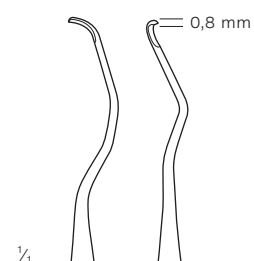

Premolars, molars, distal,
first shaft longer, working
end shorter

1/1

SP29-152-14

Fig. 13/14

17,5 cm / 6 7/8"

hollow handle, 10 mm

Ti

Gracey standard curettes, subgingival

Gracey Standard-Küretten, subgingivaler Einsatz

Gracey curetas estándar, subgingival

Gracey curettes standard, subgingival

Gracey curette standard, subgingivale


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia

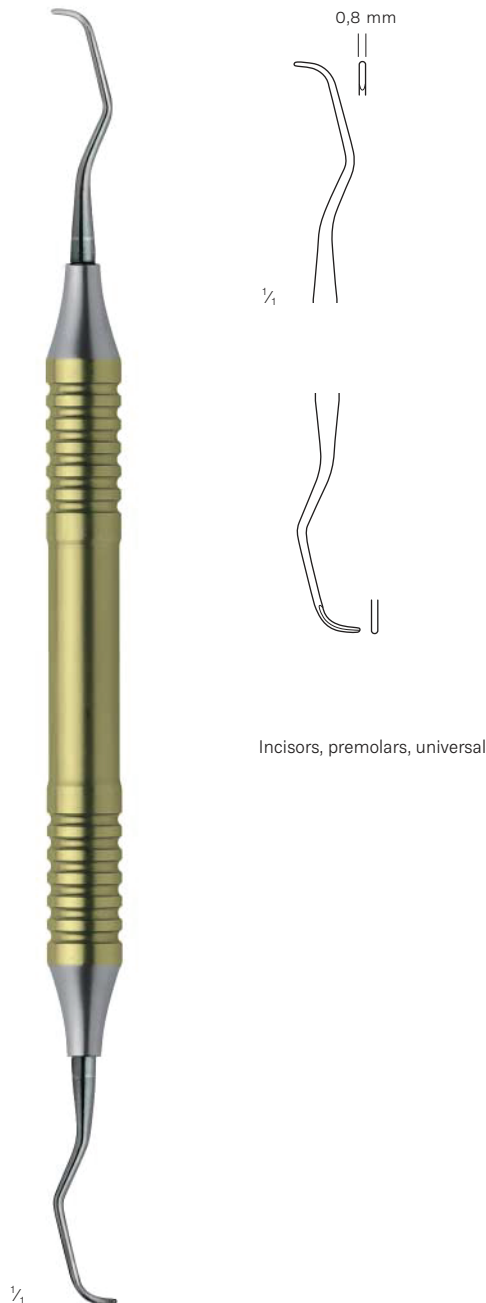


SP29-153-06

Ti

Fig. 5/6

17,5 cm / 6 7/8"

hollow handle, 10 mm

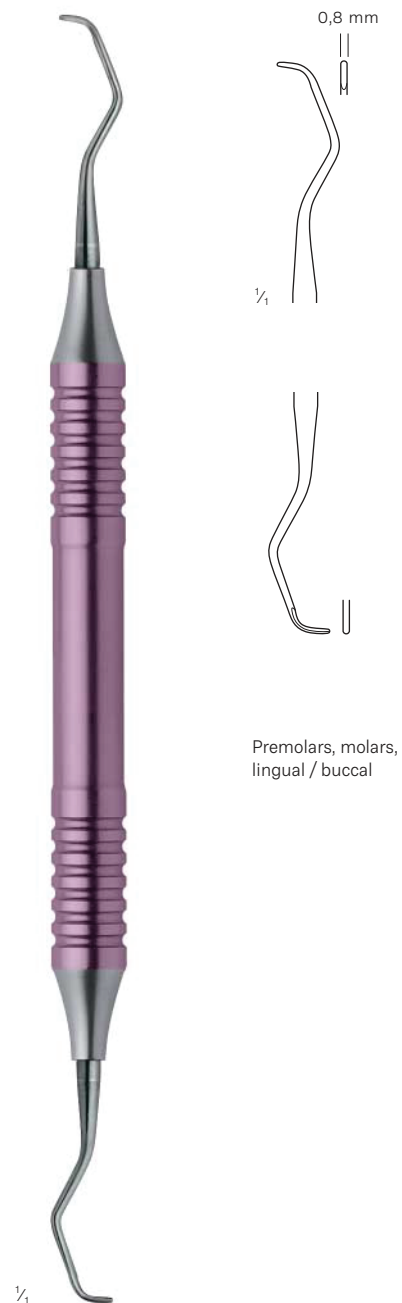


SP29-153-08

Ti

Fig. 7/8

17,5 cm / 6 7/8"

hollow handle, 10 mm

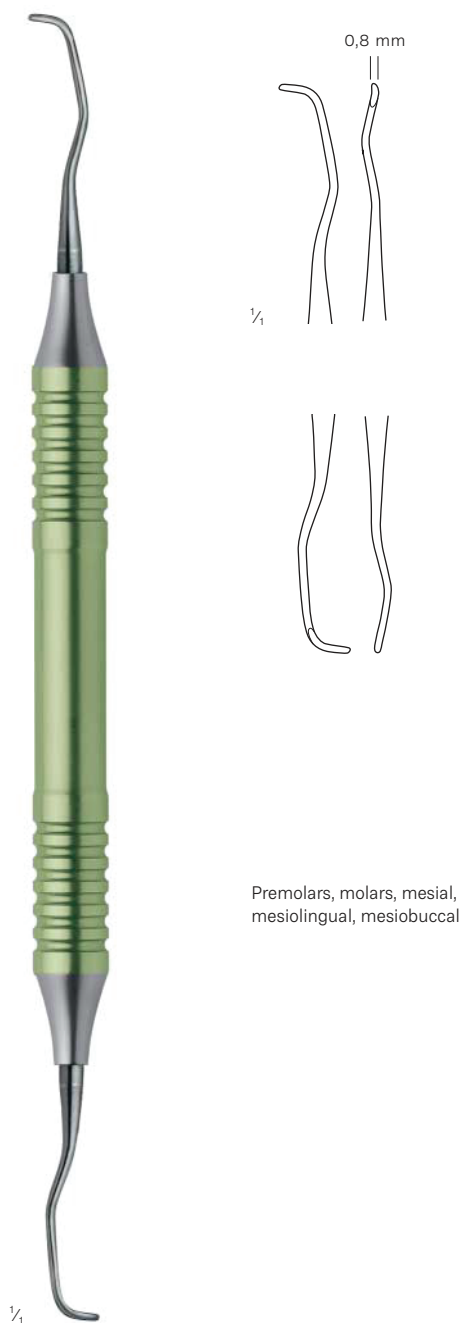


SP29-153-12

Ti

Fig. 11/12

17,5 cm / 6 $\frac{7}{8}$ "

hollow handle, 10 mm

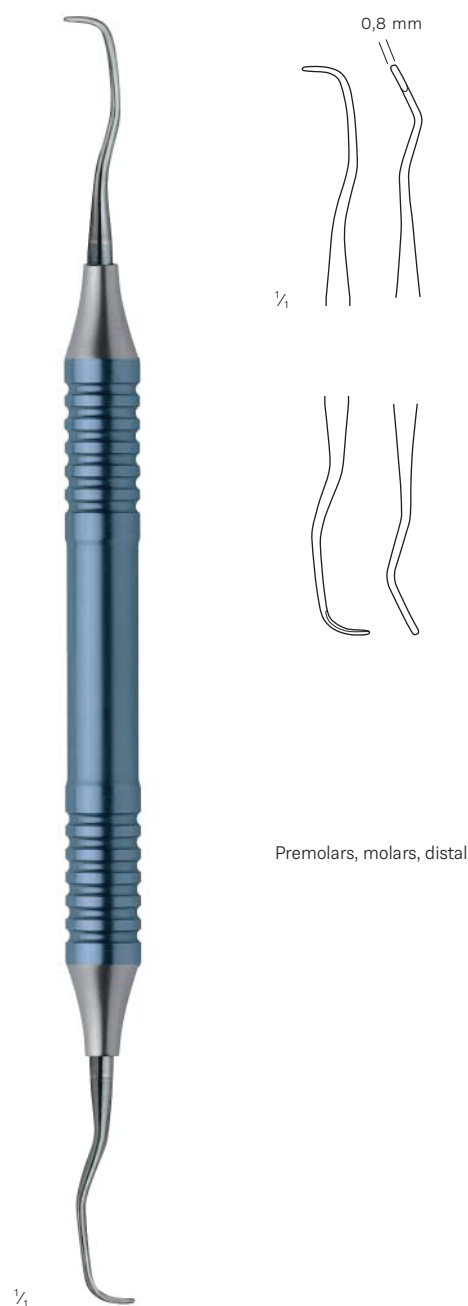

SP29-153-14

Ti

Fig. 13/14

17,5 cm / 6 $\frac{7}{8}$ "

hollow handle, 10 mm

Sharpening stones

Schleifsteine
Piedras de afilar
Meules
Pietre per affilare

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia



Record Arkansas SP29-154-00

115 x 40 x 10/2 mm
multi shape, white

For pregrinding of blunt instruments. The rounded edges and bevelled surfaces enable the sharpening of radius and surfaces



Original Arkansas SP29-154-01

100 x 25 x 7 mm
flat, grey-white

For repeated grinding and finishing of scalers and curettes



Original Arkansas SP29-154-02

100 mm lang, Ø 8 mm
round, grey-white

For burr removing from cutting edges



India-Korund SP29-154-05

115 x 45 x 2 - 10 mm

delicate, brown-beige
Para preafilar y volverles a dar forma a instrumentos



India-Korund SP29-154-06

100 mm lang, Ø 5 - 10 mm

round, delicate, brown-beige
For burr removing from cutting edges



SP29-154-08

75 mm lang, Ø 6,0 mm

Test stick plastic
package of 6 pieces
For testing the sharpness of scalers and curettes



Record Arkansas

SP29-155-00

150 x 50 mm
in wooden case



SP29-154-09

Sharpening stone oil

Periodontal probes

Parodontometer
Sondas parodontales
Sondes parodontales
Sonde parodontali

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia




SP29-156-01

Iglhaut
SP29-156-01
17,5 cm / 6 7/8"
solid handle, 6 mm



SP29-157-40

Nabers
SP29-157-40
17,5 cm / 6 7/8"
hollow handle, 8 mm
Furcation probe, graduated working end

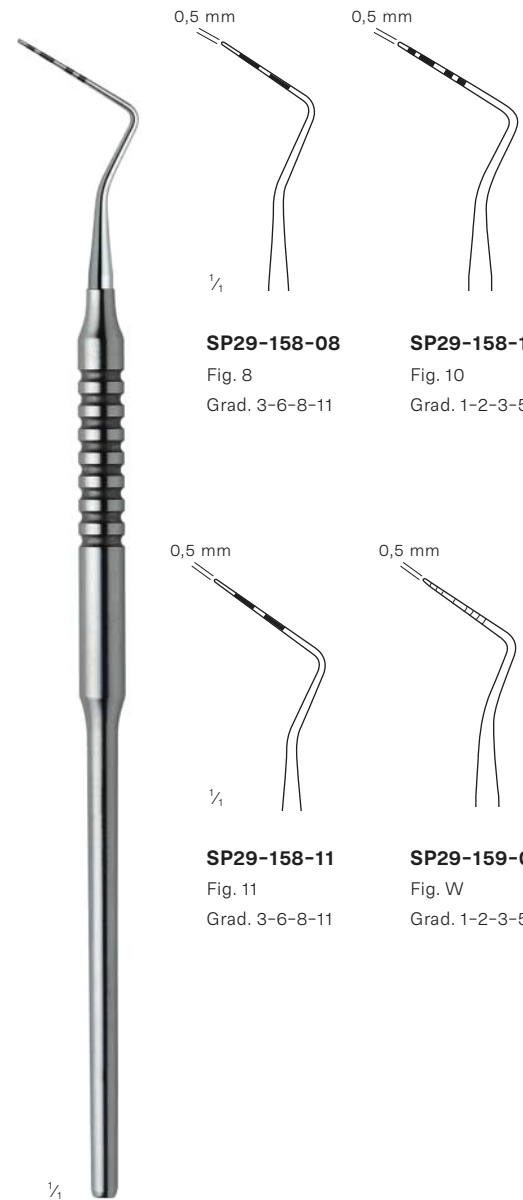

SP29-158-08

Fig. 8
Grad. 3-6-8-11

SP29-158-10

Fig. 10
Grad. 1-2-3-5-7-8-9-10

SP29-158-11


Fig. 11
Grad. 3-6-8-11

SP29-159-00

Fig. W
Grad. 1-2-3-5-7-8-9-10

SP29-158-08 - SP29-159-00

16 cm / 6 1/8"
solid handle, 6 mm

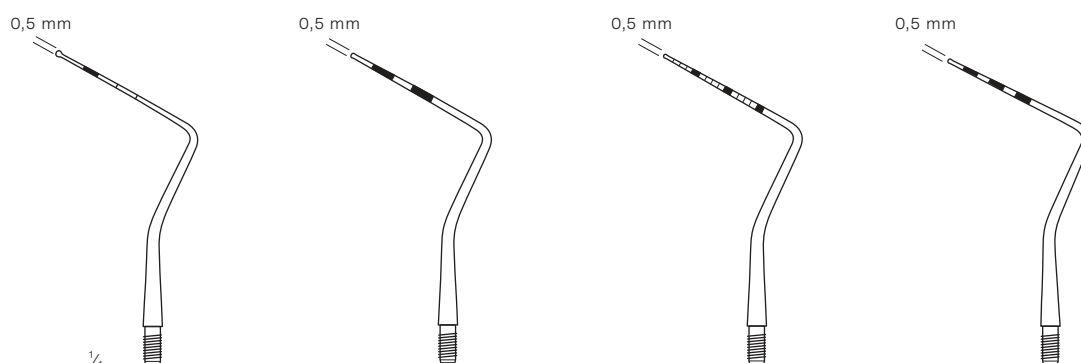

SP29-160-11	SP29-160-13	SP29-160-21	SP29-160-23
SP29-161-11	SP29-161-13	SP29-161-21	SP29-161-23
Fig. CPG 11,5 WHO Grad. 3,5-5,5-8,5-11,5 with ball	Fig. CPG 12 Grad. 3-6-9-12	Fig. CP 15 UNC Grad. 1-2-3-4-...-14-15	Fig. CPNG 22 Grad. 2-4-6-8-10-12




Probe handles see page 244

Periodontal probes

Parodontometer
Sondas parodontales
Sondes parodontales
Sonde parodontali

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia

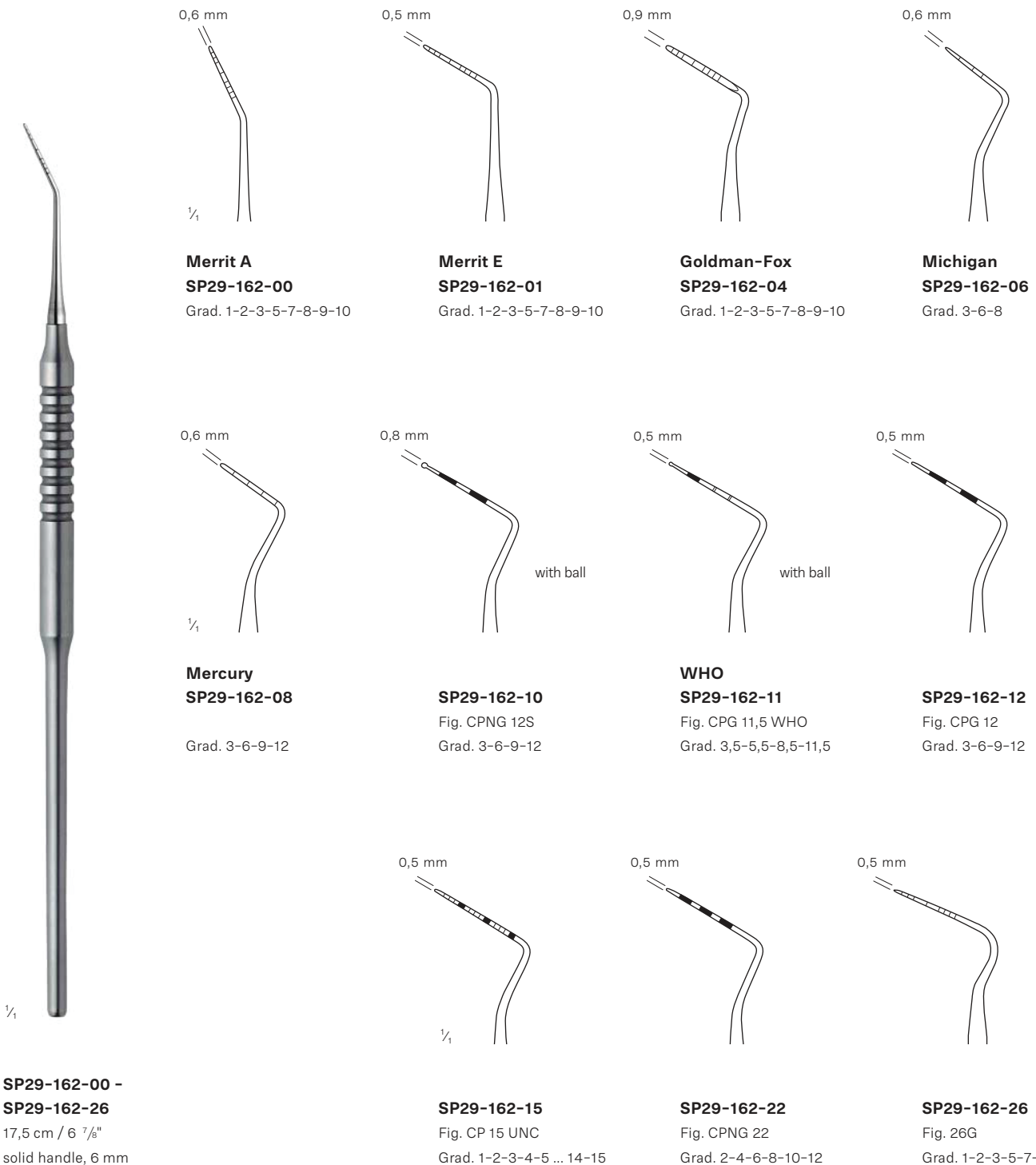



0,8 mm

1/1

Williams-Goldman-Fox

SP29-163-02

Grad. 1-2-3-5-7-8-9-10

0,6 mm

SP29-163-04

Fig. CP 11/23

Grad. 3-6-8-11

0,6 mm

SP29-163-06

Fig. CP 12/23

Grad. 3-6-9-12

0,5 mm

0,8 mm

SP29-163-08

Fig. GF/3

Grad. 1-2-3-5-7-8-9-10





0,5 mm

0,8 mm

SP29-163-10

Fig. UM/3

Grad. 3-6-8

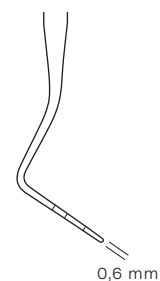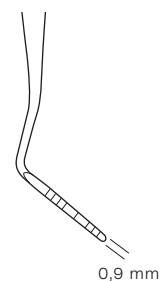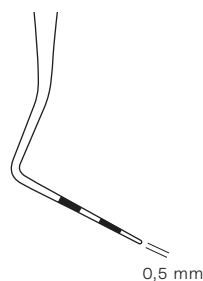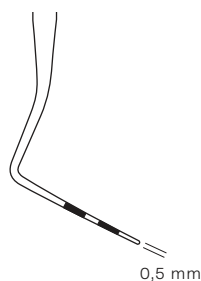

0,5 mm

0,8 mm

1/1

SP29-163-12

Fig. W/3

Grad. 1-2-3-5-7-8-9-10

0,6 mm

SP29-163-14

Fig. GF/6

Grad. 1-2-3-5-7-8-9-10

0,6 mm

SP29-163-16

Fig. UM/6

Grad. 3-6-8

0,5 mm

SP29-163-18

Fig. W/6


Grad. 1-2-3-5-7-8-9-10

0,6 mm


SP29-163-20

Fig. GF/17


Grad. 1-2-3-5-7-8-9-10

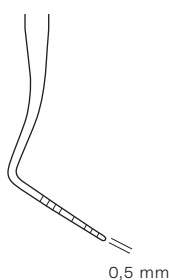


0,5 mm

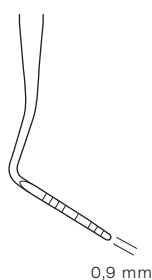


0,9 mm

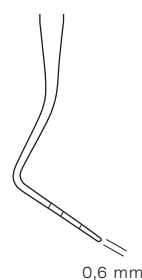

0,6 mm

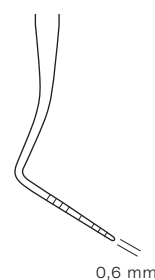

0,6 mm

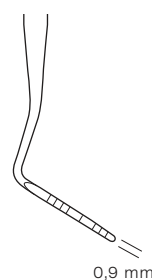

0,9 mm

**SP29-163-02 -
SP29-163-20**
17,5 cm / 6 7/8"
hollow handle, 8 mm

Periodontal probes

Parodontometer
Sondas parodontales
Sondes parodontales
Sonde parodontali

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia



0,5 mm



1/4

SP29-163-22

Fig. UM/17

Grad. 3-6-8



0,5 mm

0,5 mm



SP29-163-24

Fig. W/17

Grad. 1-2-3-5-7-8-9-10



0,5 mm

0,9 mm



SP29-163-26

Fig. GF/23

Grad. 1-2-3-5-7-8-9-10



0,6 mm

0,6 mm



SP29-163-28

Fig. UM/23

Grad. 3-6-8



0,5 mm

0,6 mm



SP29-163-30

Fig. W/23

Grad. 1-2-3-5-7-8-9-10



0,5 mm


SP29-163-22 - SP29-163-30

17,5 cm / 6 7/8"

hollow handle, 8 mm




Cross
SP29-164-10
17,5 cm / 6 7/8"
hollow handle, 8 mm

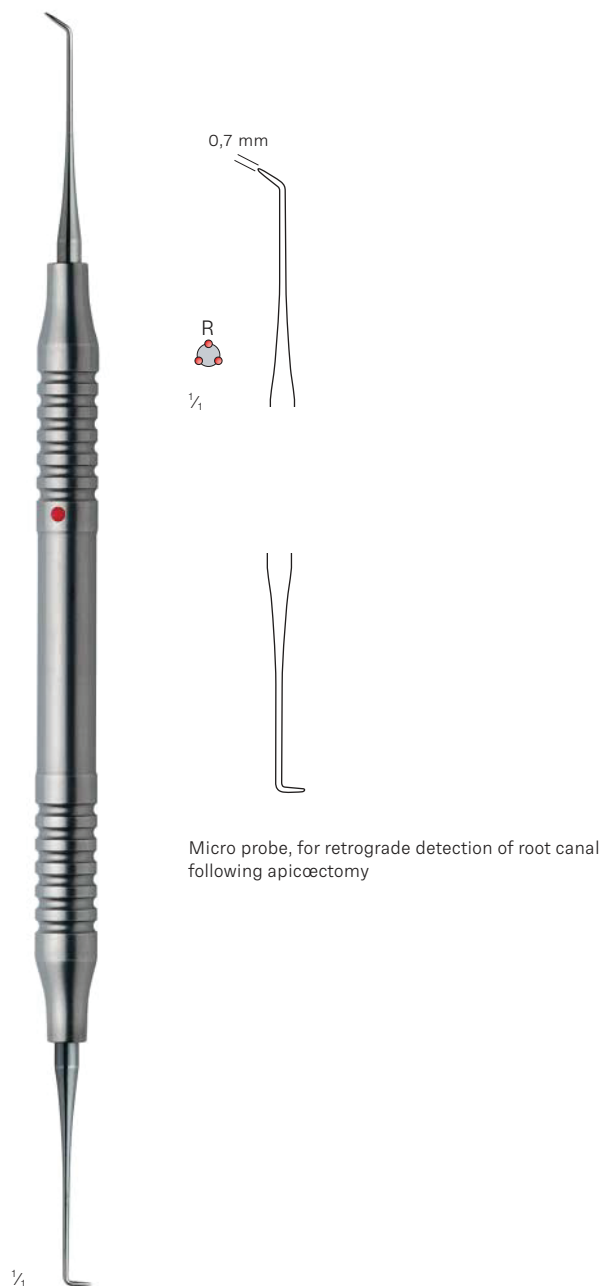

Zürich
SP29-165-02
17,5 cm / 6 7/8"
hollow handle, 8 mm




FinoLine

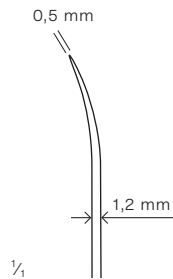

SP29-165-03
Fig. 3




SP29-165-06
Fig. 6

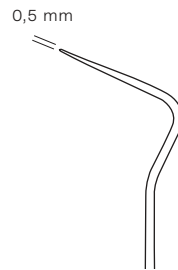

SP29-165-08
Fig. 8 A




SP29-165-09
Fig. 9

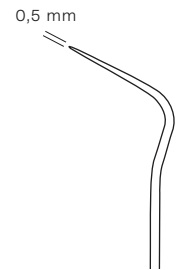

SP29-165-12
Fig. 12




SP29-165-17
Fig. 17

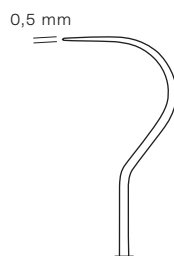


SP29-165-23
Fig. 23



SP29-165-30
Fig. 3A

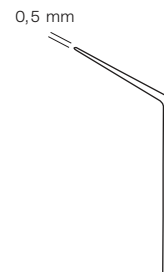


SP29-165-60
Fig. 6XL

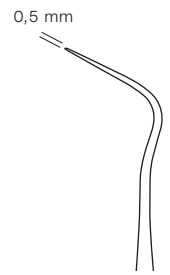

White
SP29-165-12w
Fig. 12

SP29-165-03 -
SP29-165-12w
17,5 cm / 6 7/8"
extra delicate tooth probe
solid handle, 6 mm

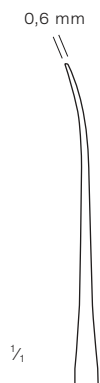

SP29-166-03

Fig. 3


Detection of furcations and for easy entry into interproximal spaces



SP29-166-06


Fig. 6

Detection of plaque, caries, root canals and furcations

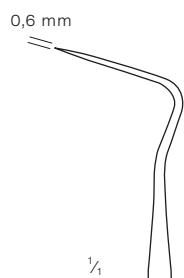

SP29-166-08

Fig. 8A

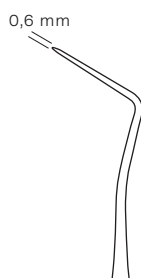


SP29-166-09

Fig. 9



Briault
SP29-166-11

Fig. EXD 11

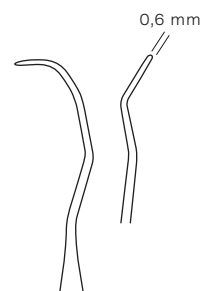


SP29-166-12

Fig. EXD 11-12


SP29-166-03 -
SP29-166-12

17,5 cm / 6 7/8"
solid handle, 6 mm

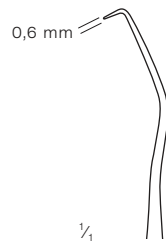


White
SP29-166-17
Fig. 17

Mainly for sub-gingival examination and approx. caries

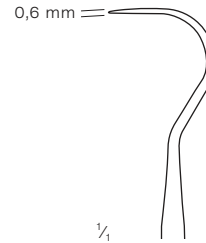


White
SP29-166-23
Fig. 23


For general examination, especially for incisors

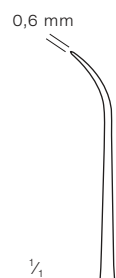

SP29-166-33
Fig. 3A

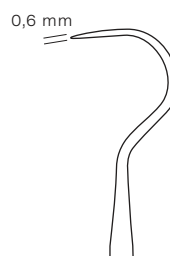

SP29-166-54
Fig. 54

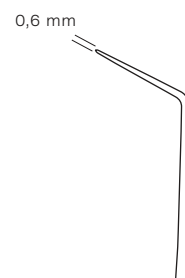


SP29-166-66
Fig. 6XL

Detection of plaque, caries, root canals and furcations

SP29-166-17 -
SP29-166-66
17,5 cm / 6 7/8"
solid handle, 6 mm





**SP29-167-01 -
SP29-167-23**
17,5 cm / 6 7/8"
hollow handle, 8 mm

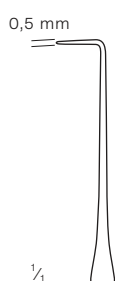

SP29-167-01
Fig. UT 1/2



SP29-167-02
Fig. UT 2/2

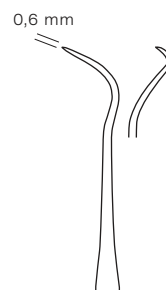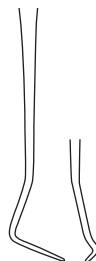



SP29-167-03
Fig. 3CH/3CH





SP29-167-04
Fig. 6/3




SP29-167-06
Fig. 23/6

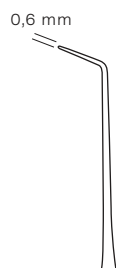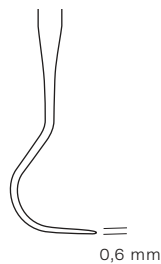


SP29-167-16
Fig. 16/17

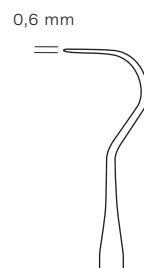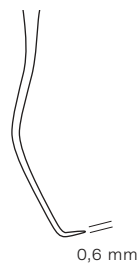


SP29-167-17
Fig. 17/23

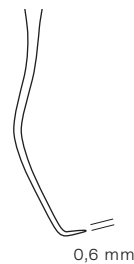


SP29-167-23
Fig. 23/17A

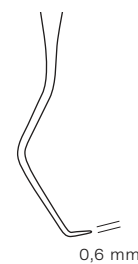

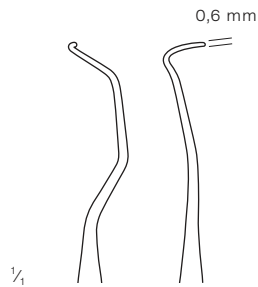


Nabers
SP29-168-01
Fig. 1

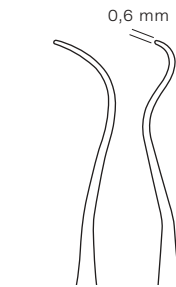





Nabers
SP29-168-02
Fig. 2




SP29-168-05
Fig. 5/6

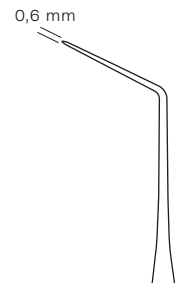

SP29-168-23
Fig. 16/23

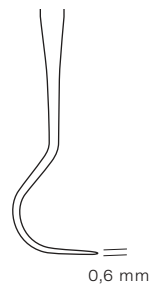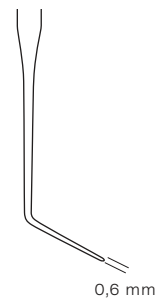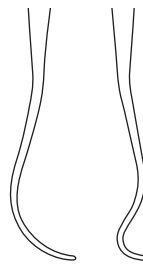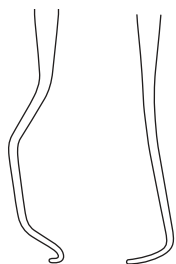

SP29-168-01 -
SP29-168-23
17,5 cm / 6 7/8"
hollow handle, 8 mm

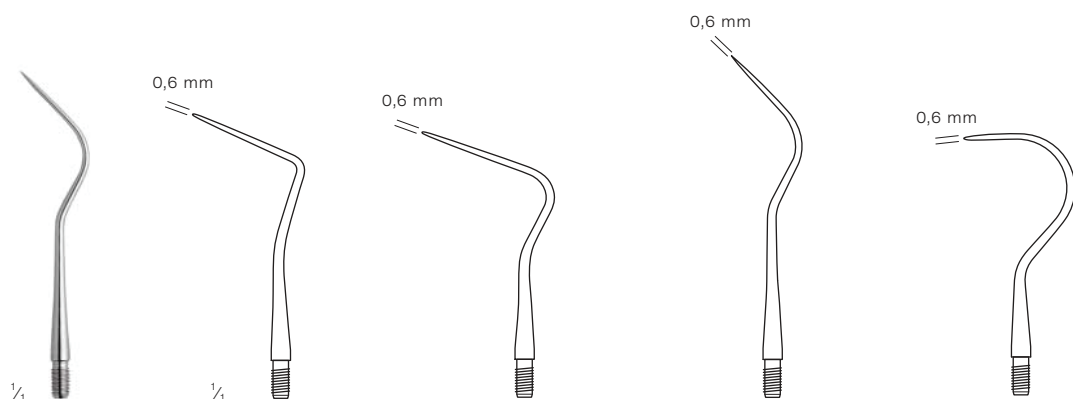

SP29-169-08

Fig. EXG 8

SP29-169-09

Fig. EXG 9

SP29-169-10

Fig. EXG 10

SP29-169-23

Fig. EXG 23




1/4


Probe handles

Excavators
 Exkavatoren
 Excavadores
 Excavateurs
 Escavatori

Scalers, periodontology
 Reinigungsinstrumente, Parodontologie
 Limpiadores de sarro, parodontología
 Instruments à détartrer, parodontologie
 Ablatori di tartaro, parodontologia

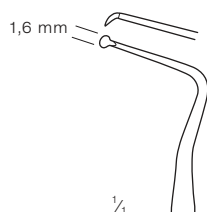


SP29-170-01
 Fig. GM 1A

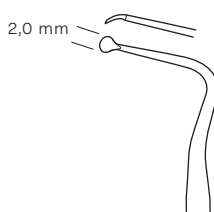


SP29-170-02
 Fig. GM 2A



SP29-171-00
 Fig. GM 0

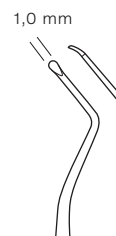

SP29-171-01
 Fig. GM 1

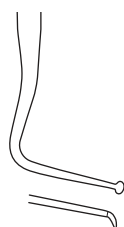

**SP29-170-01 -
 SP29-171-34**
 17,5 cm / 6 7/8"
 hollow handle, 8 mm
 double ended



grey



yellow



red



green



violet




white



blue




light-blue

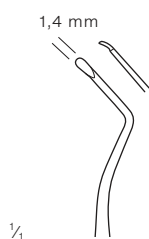

SP29-171-02

Fig. GM 2

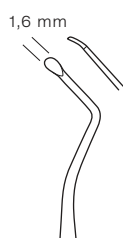


SP29-171-03

Fig. GM 3



SP29-171-04

Fig. GM 4

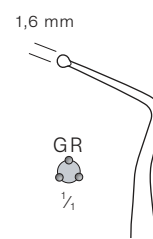

Zürich

SP29-171-34

Ø 1,9 mm





1/1

**SP29-172-00 -
 SP29-172-44**
 17,5 cm / 6 7/8"
 hollow handle, 8 mm
 double ended




SP29-172-00
 Fig. GM 125/126



SP29-172-02
 Fig. GM 127/128




SP29-172-04
 Fig. GM 129/130

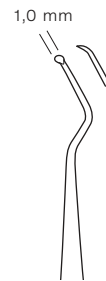

SP29-172-06
 Fig. GM 133/134



SP29-172-08
 Fig. GM 153/154

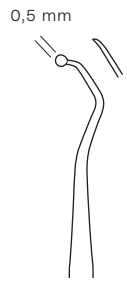

SP29-172-10
 Fig. GM 14




SP29-172-12
 Fig. GM 38/39




SP29-172-14
 Fig. GM E1



SP29-172-16
 Fig. GM E2

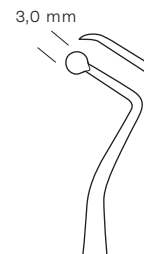

SP29-172-18
 Fig. GM E3

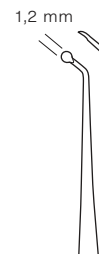

SP29-172-20
 Fig. GM 6



SP29-172-22
 Fig. GM 17S






SP29-172-24
Fig. GM 17



SP29-172-26
Fig. GM 18

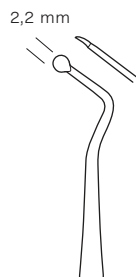

SP29-172-28
Fig. GM 19



SP29-172-30
Fig. GM 17W



SP29-172-32
Fig. GM 18W



SP29-172-34
Fig. GM 19W




SP29-172-35
Fig. GM 61/62




SP29-172-36
Fig. GM 63/64




SP29-172-38
Fig. GM 65/66

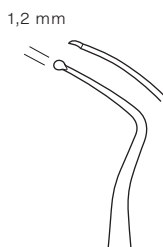


SP29-172-40
Fig. GM 31L

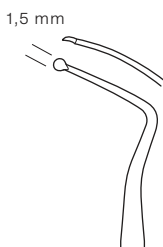

SP29-172-42
Fig. GM 32L

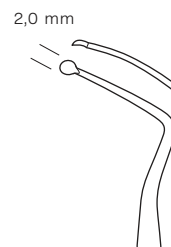

SP29-172-44
Fig. GM 33L

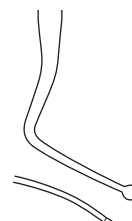

Enamel chisels

Schmelzmeißel
Cinceles para esmalte
Ciseaux à émail
Tagliasmalto

Scalers, periodontology

Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia



1/1

Black


**SP29-173-01 -
SP29-173-80**

17,5 cm / 6 7/8"
solid handle, 6 mm
double ended




SP29-173-01

Fig. 1

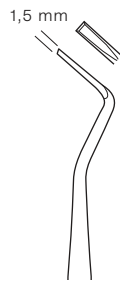


SP29-173-02

Fig. 2

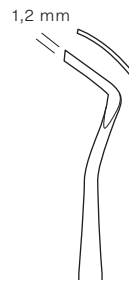

SP29-173-26

Fig. 26

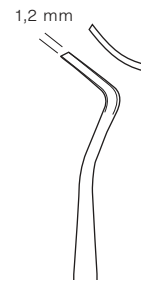

SP29-173-27

Fig. 27



SP29-173-28

Fig. 28



SP29-173-29

Fig. 29



SP29-173-52

Fig. 51/52



SP29-173-54

Fig. 53/54



SP29-173-64

Fig. 63/64



SP29-173-66

Fig. 65/66



SP29-173-78

Fig. 77/78




SP29-173-80

Fig. 79/80




1/1

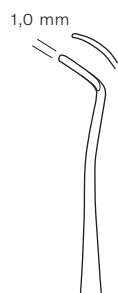






Wedelstaedt
SP29-174-02 -
SP29-174-06
17,5 cm / 6 7/8"

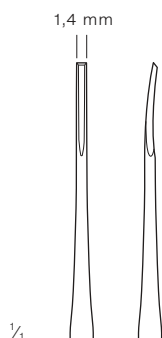


SP29-174-02
Fig. 1/2

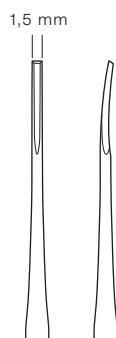

SP29-174-04
Fig. 3/4

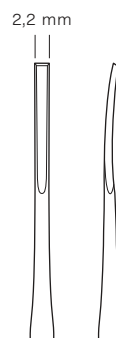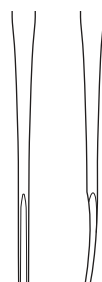

SP29-174-06
Fig. 5/6

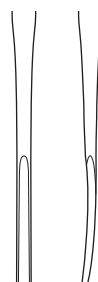

solid handle, 6 mm
double ended

Root canal instruments

Nervkanal-Instrumente

Instrumentos para canales

Instrumenti a canaux

Strumenti per trattamenti canalari


Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instrumenti a détartrer, parodontologie

Ablatori di tartaro, parodontologia

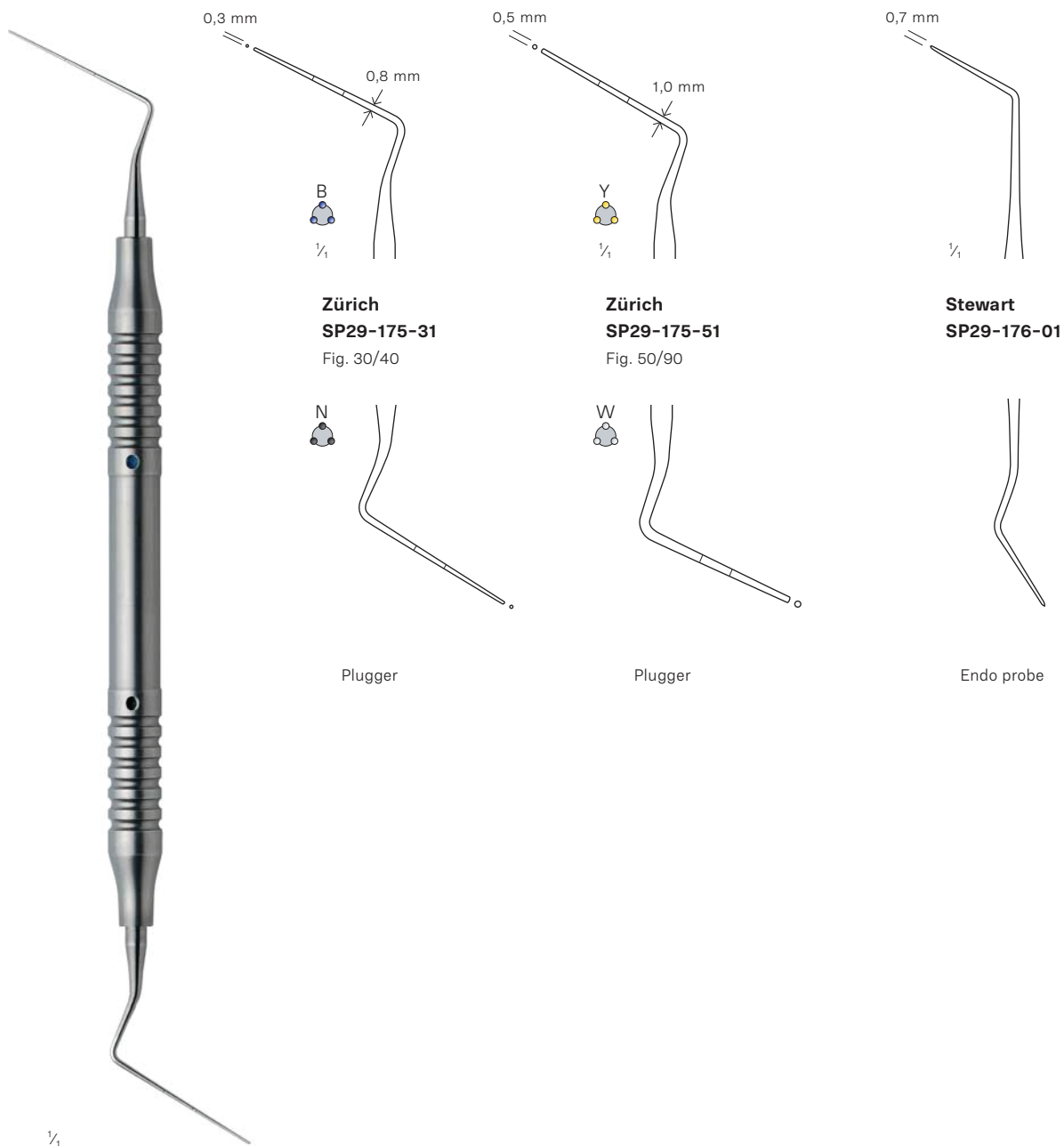

SP29-175-31 -


SP29-176-01

17,5 cm / 6 7/8"

hollow handle, 8 mm


double ended

GR grey Y yellow R red G green V violet W white B blue N black

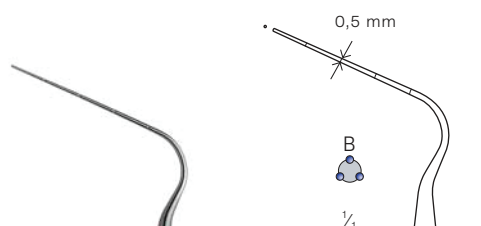


SP29-177-03

Fig. 30

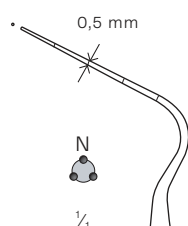


SP29-177-04

Fig. 40

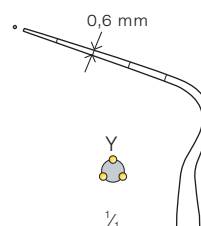


SP29-177-05

Fig. 50

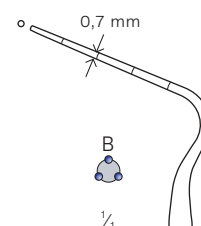


SP29-177-06

Fig. 60

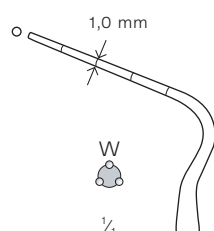


SP29-177-09

Fig. 90

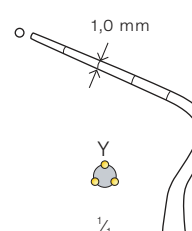

SP29-177-10

Fig. 100

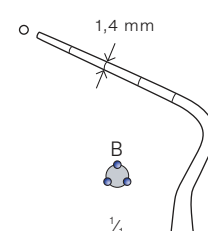

SP29-177-12

Fig. 120

Plugger



1/4

SP29-177-03 -

SP29-177-12


17,5 cm / 6 7/8"
solid handle, 6 mm
single ended

Root canal instruments


Nervkanal-Instrumente
Instrumentos para canales
Instruments à canaux
Strumenti per trattamenti canalari

Scalers, periodontology


Reinigungsinstrumente, Parodontologie
Limpiadores de sarro, parodontología
Instruments à détartrer, parodontologie
Ablatori di tartaro, parodontologia

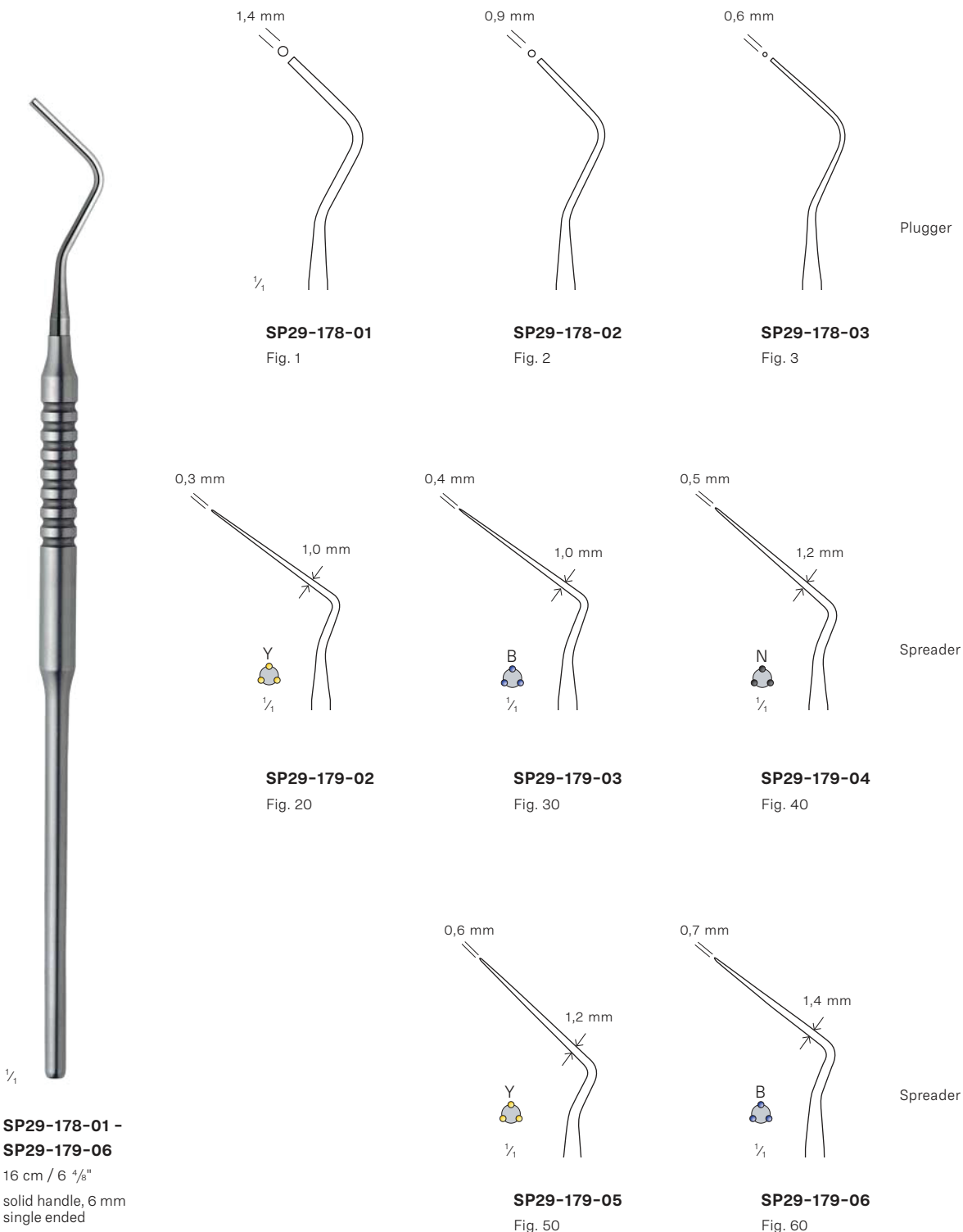

GR grey Y yellow R red G green V violet W white B blue N black

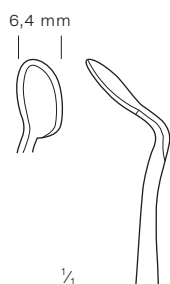

SP29-180-37
Fig. 7

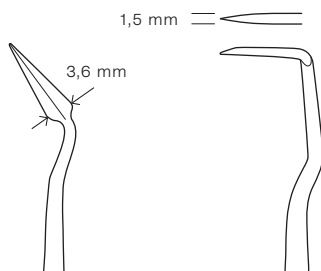

SP29-180-38
Fig. 8



SP29-180-39
Fig. 9



SP29-180-40
Fig. 10



SP29-180-41
Fig. 11



**SP29-180-37 -
 SP29-180-41**
 17,5 cm / 6 7/8"
 solid handle, 6 mm
 double ended

Gingivectomy knives

Gingivektomiemesser

Bisturís para la gingivectomía

Bistouris pour la gingivectomie

Bisturie per gengivectomia

Scalers, periodontology

Reinigungsinstrumente, Parodontologie

Limpiadores de sarro, parodontología

Instruments à détartrer, parodontologie

Ablatori di tartaro, parodontologia



1/4

Allen-Iglhaut


SP29-181-03 -

SP29-181-04



16 cm / 6 4/8"

hollow handle, 10 mm


double ended

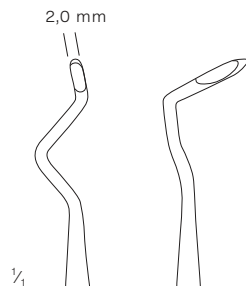


SP29-181-03

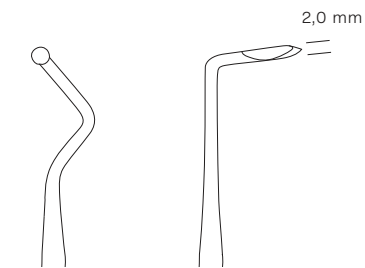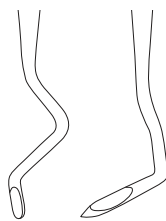


SP29-181-04

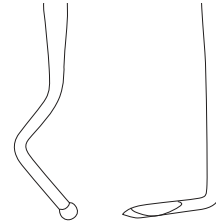

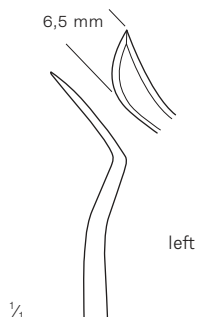


Kirkland
SP29-182-51


For gingivectomy or gingivoplastic surge
 Very effective in retro molar region

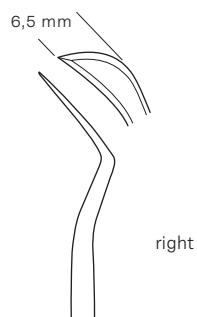

Kirkland
SP29-182-61

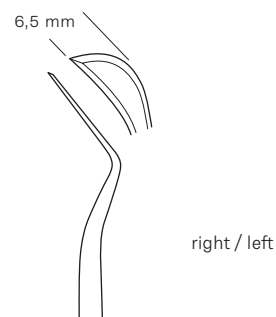


Kirkland
SP29-183-61

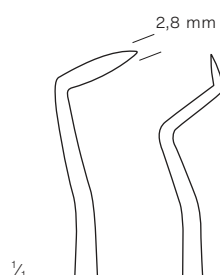


Orban
SP29-184-01

Fig. 1

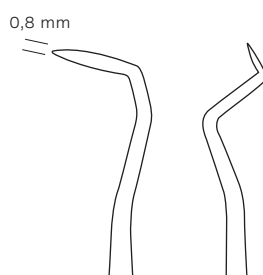

Orban
SP29-184-02

Fig. 2

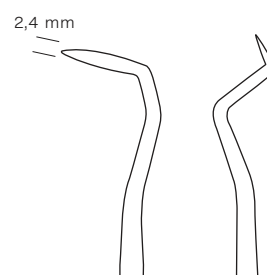

Orban
SP29-184-12

Fig. 1/2



SP29-182-51 -
SP29-184-12
 16 cm / 6 1/8"
 solid handle, 6 mm
 single ended

Counter angled shank, especially suitable for excision of
 interproximal tissue in gingivectomy